

GENIUS

Pure precise
coagulation
every time.

CYNOSURE®
LUTRONIC®






AI meets RFMN

Only Genius uses an AI-based Radio Frequency MicroNeedling (RFMN) system that ensures **pure and precise, deep coagulation** every time.

Genius utilizes **real time tissue feedback** to continually adjust energy and create thousands of precise & consistent microcoagulation zones.

This effectively stimulates new collagen and elastin for consistent volumization. **Genius is safe for all skin types, all year round.**

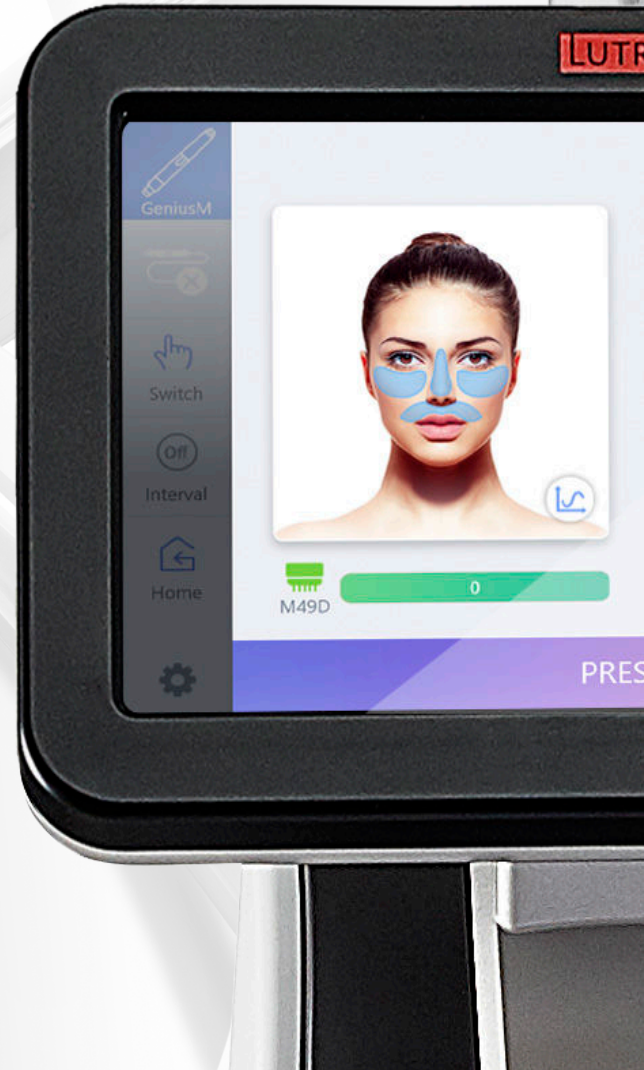



Real time tissue sensing and energy adjustment /

- The only RFMN system to monitor and adjust energy in real time, over **500 times per second**
- Precise depth of penetration for safe treatments and **clinical confidence** in all patient skin types
- **Efficacious face/body volumization** for the appearance of tighter skin
- Measured energy delivery for **higher quality results** and ROI
- Integrated air-cooling capability available
- Complementary **treatment versatility**. Ultra Skin Solutions™



Real patients. Genuine concerns/

SKIN REVITALIZATION Courtesy of Omar Ibrahim, MD, PHD



Before and after 3 treatments

SKIN REVITALIZATION Courtesy of Steven Weiner, MD



Before and after 3 treatments

SKIN REVITALIZATION Courtesy of Jason Pozner, MD



Before and after 2 treatments

SKIN REVITALIZATION Courtesy of Lasting Beauty Medspa



Before and after 2 treatments

Life-changing results/

SKIN REVITALIZATION Courtesy of Jason Pozner, MD



Before and after 2 treatments

SKIN REVITALIZATION Courtesy of A. Jay Burns, MD




Before and after 2 treatments

SKIN REVITALIZATION Courtesy of Paul M. Friedman, MD



Before and after 3 treatments

SKIN REVITALIZATION Courtesy of Swapnil Shah, MD

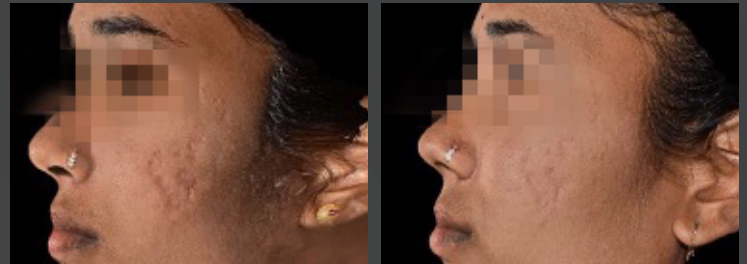


Before and after 3 treatments

Confidence


RFMN is a blind treatment. Without advanced tissue feedback technology, other RFMN devices cannot know what is happening under the skin.

Competitor 1

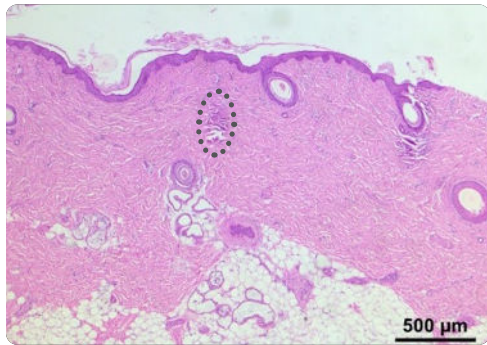


Minimal Tissue injury from uncontrolled low energy delivery

Competitor 2

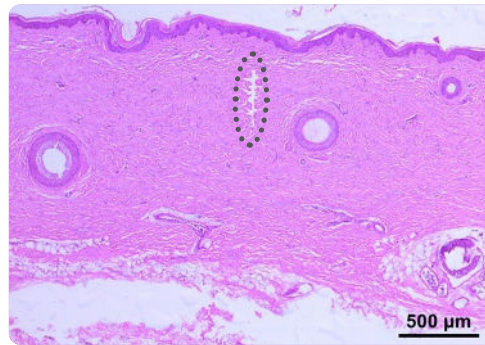

Tissue injury from uncontrolled high power = Desiccation

Genius

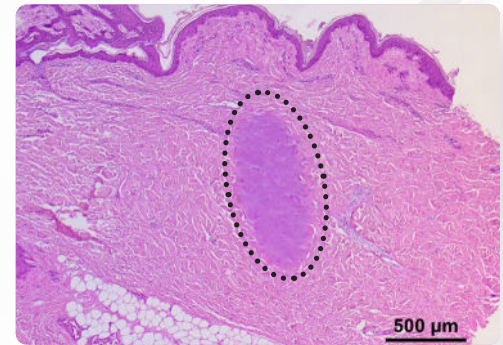

Full and Exact Energy Delivered
Pure Coagulation

We're your trusted *partner on this journey/*

Technology

- Forward thinking science
- Advanced AI-based RF Microneedling

Product

- Significant clinical evidence
- Highest quality technological improvements
- Only company to manufacture its own sharper, thinner needles

Clinical

- Comprehensive virtual training
Personalized in-office training for
beginner and advanced teams

Service

- Experienced technical support
team 24/7
- Expert field service team

Practice Development

- Industry expert consultants
- Focused on growing your business

Amplified Marketing Practice Support

- Dedicated patient support team
- Comprehensive turnkey marketing





Visit/
www.lutronic.com

Follow us/


@cynosurelutronicemea Cynosure Lutronic EMEA

Intended for distribution in US only.

© 2024 CYNOSURE LUTRONIC, INC. All Rights Reserved. Lutronic Genius is a trademark of Lutronic Aesthetic, Inc. Cynosure, LLC and Lutronic Aesthetic, Inc. do not diagnose or treat medical conditions, prescribe treatment protocols, provide medical advice, or otherwise practice medicine.

The Genius aesthetic workstation is indicated for the following: Dermatologic and general surgical procedures for electrocoagulation and hemostasis. CUS-1623-EMEA

CYNOSURE[™]
LUTRONIC[™]

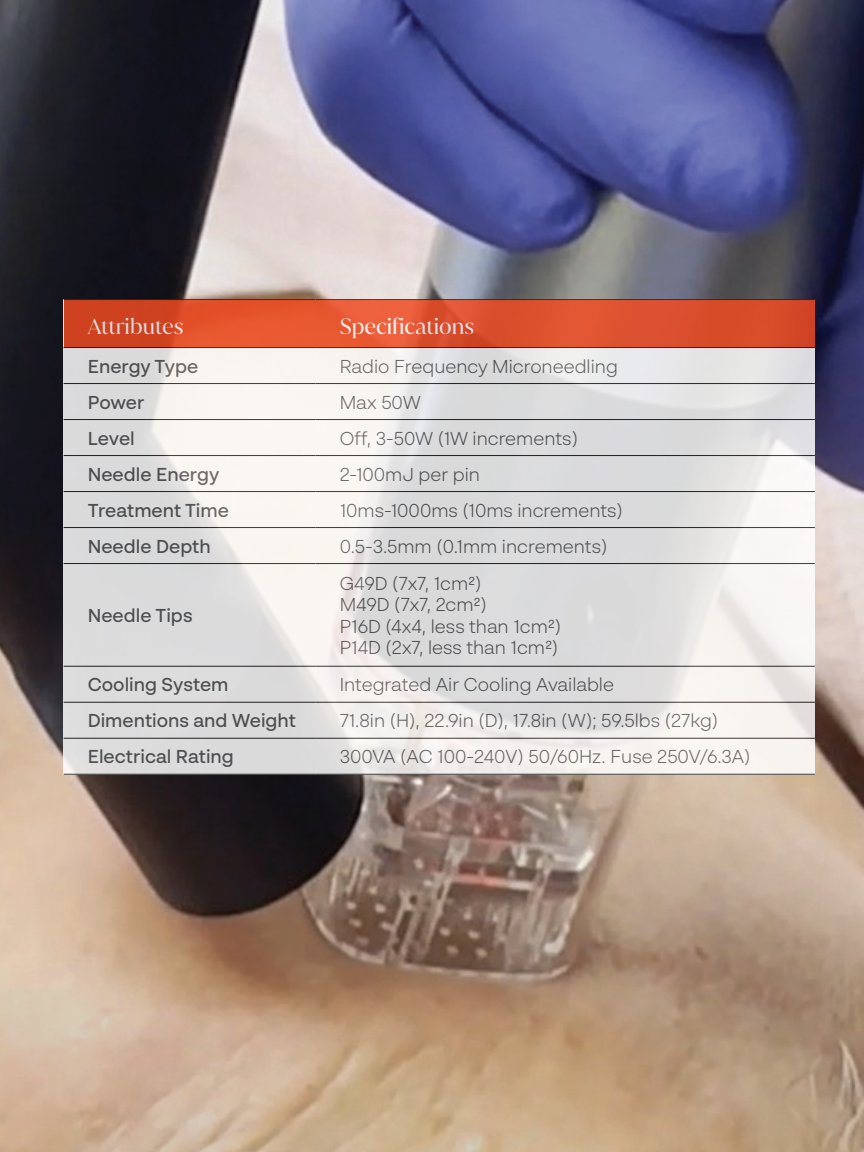

Attributes	Specifications
Energy Type	Radio Frequency Microneedling
Power	Max 50W
Level	Off, 3-50W (1W increments)
Needle Energy	2-100mJ per pin
Treatment Time	10ms-1000ms (10ms increments)
Needle Depth	0.5-3.5mm (0.1mm increments)
Needle Tips	G49D (7x7, 1cm ²) M49D (7x7, 2cm ²) P16D (4x4, less than 1cm ²) P14D (2x7, less than 1cm ²)
Cooling System	Integrated Air Cooling Available
Dimensions and Weight	71.8in (H), 22.9in (D), 17.8in (W); 59.5lbs (27kg)
Electrical Rating	300VA (AC 100-240V) 50/60Hz. Fuse 250V/6.3A